


Ligating Clips Appliers with H·E·R·O™ Pressure Control System Highly focused on patient's safety

H·E·R·O™

Patent no. 11,617,595

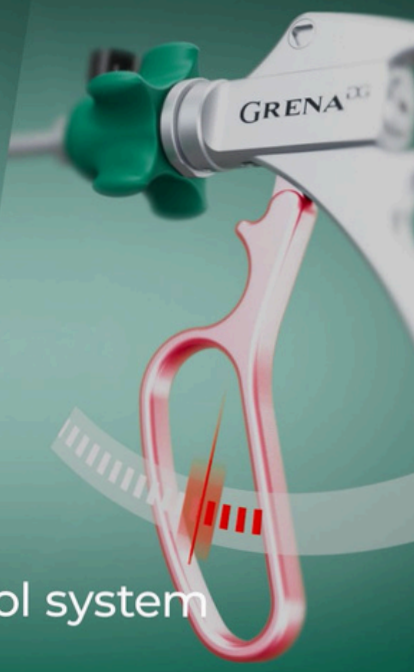

Innovative pressure control system

PROVIDING SPECIALISTS WITH
INNOVATIVE
HEALTHCARE SOLUTIONS

WITH THE CENTRAL FOCUS ON
PATIENTS' SAFETY AND USERS' COMFORT



Ligating Clips Appliers with H·E·R·O™ Pressure Control System

For more security we equipped our laparoscopic appliers with a breakthrough safety system - H·E·R·O™. H·E·R·O™ system prevents from too much force to be transferred to the applicator jaws protecting the tissues as well as the device. H·E·R·O™ stands for High Energy Override.

In some cases, too much pressure can damage blood vessel doing harm to the patient. In response we designed an innovative pressure con-

trol system so that operator can focus on applying the clip without considering the closing force.

Our new technology is a safety system that reduces the mechanical force protecting the tissue during the procedure.



H·E·R·O™

H·E·R·O™ technology protects the tissue, the clip and the applicator from damage making the procedure more secure for the patient as well as the operator.

REF	CLIP SIZE	TYPE	WORKING LENGTH (CM)	SHAFT/TROCAR DIAMETER (MM)	COLOUR	PCS PER BOX
0301-04MEN	M	Endoscopic	33	5/5	●	1
0301-04MEBN	M	Endoscopic	45	5/5	●	1
0301-04MLEN	ML	Endoscopic	33	5/5	●	1
0301-04MLEBN	ML	Endoscopic	45	5/5	●	1
0301-04LEN	L	Endoscopic	33	10/10	●	1
0301-04LEBN	L	Endoscopic	45	10/10	●	1
0301-04XLEN	XL	Endoscopic	33	10/10	●	1
0301-04XLEBN	XL	Endoscopic	45	10/10	●	1
0301-04XXLEN	XXL	Endoscopic	33	10/10	●	1
0301-04XXLEBN	XXL	Endoscopic	45	10/10	●	1

Contact Austria:

HCP Healthcarepartner GmbH
Wolfgangsee Straße 26, 5023 Salzburg-Guggenthal, Österreich
T +43 (0)662/62 12 43 • F -11 • office@hcp-austria.com

hcp-austria.com